

MVE HEco™ 800 Series Freezers




Features include

- LN₂ efficiency
- Fully enclosed wiring and plumbing
- Vacuum jacketed transfer hose
- Dry sample storage
- -190°C top box temperature
- Lowest liftover height
- Single step standard
- Optional double-tier step



- Standard One-Tier Step
- Optional Two-Tier Steps

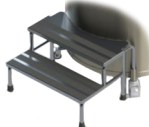

Dura-Cyl® 230
Liquid Cylinders
for LN₂ use




Featuring the Chart MVE Touch Screen & TEC 3000

The MVE HEco Series is the next generation High Efficiency freezer from Chart MVE. With their unique shroud design, streamlined LN₂ plumbing, and vacuum jacketed transfer hose, the MVE HEco 800 Series freezers provide efficient use of LN₂, making them the most efficient vapor freezers available. These next generation high efficiency freezers incorporate hinged work surfaces that fully enclose all electronics and plumbing to enhance overall safety and usability.

The technical and visual improvements combine to create a more aesthetically pleasing freezer that offers greater functionality. As with all of MVE's cryogenic freezers, the MVE HEco 800 Series freezers are available in three unique sizes and provide maximum storage density as well as the longest hold time and highest sample security in the industry.

MVE HEco 800 Temperature Test*

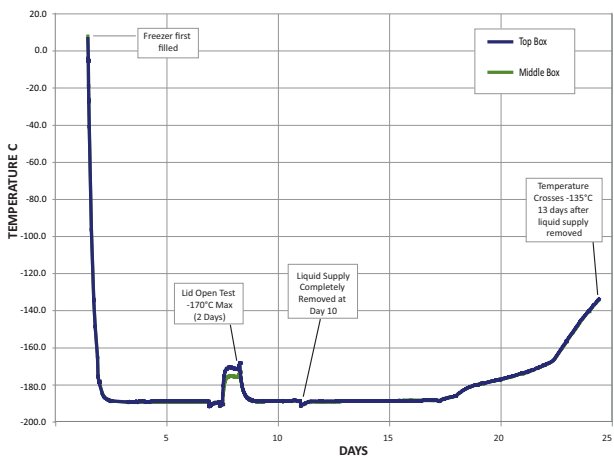

Temperature Test Graph

* Temperature test indicates typical performance of HEco Series freezer with full inventory system and factory recommended level settings. Actual performance may vary with atmospheric conditions and usage.

www.chartbiomed.com
844-MVE-CRYO



	MVE HEco 815P-190	MVE HEco 818P-190	MVE HEco 819P-190						
Maximum Storage Capacity									
1.2 & 2 ml Vials (Internally Threaded)	15,600	18,200	19,500						
Quantity of Large Racks 100 cell boxes	12	12	12						
Quantity of Mini Racks 25 cell boxes	4	4	4						
Number of Shelves per Rack	12	14	15						
Performance									
LN ₂ Capacity w/o Inventory L est.	349	403	446						
LN ₂ Capacity at Vapor Platform L est.	51	55	55						
Unit Dimensions									
Neck Opening in. (mm)	12.5 (317)	12.5 (317)	12.5 (317)						
Usable Internal Height in. (mm)	26.5 (673)	30.7 (781)	34.5 (877)						
Inner Diameter in. (mm)	28.8 (731)	28.8 (731)	28.8 (731)						
Overall Height in. (mm)	49.5 (1257)	53.9 (1370)	57.70 (1465)						
Liftover Height in. (mm)	39 (991)	43.8 (1115)	47.63 (1209)						
Door Width Requirement** in. (mm)	32.0 (813)	32.0 (813)	32.0 (813)						
Door Width Requirement, with handles** in. (mm)	33.2 (843)	33.2 (843)	33.2 (843)						
Weight Empty* lb. (kg)	475 (215)	495 (225)	515 (234)						
Weight Liquid Full* lb. (kg) est.	1134 (514)	1168 (530)	1340 (608)						
Blood Bag Capacities									
	Total Bags	Bags/ Frame	No. Frames	Total Bags	Bags/ Frame	No. Frames	Total Bags	Bags/Frame	No. Frames
791 OS/U Medsep (25 ml)	1,224	6	204	1428	7	204	1,836	9	204
4R9951 (50 ml)	768	6	128	896	7	128	1024	8	128
4R9953 (250 ml)	416	4	104	416	4	104	520	5	104
4R9955 (500 ml)	304	4	76	304	4	76	380	5	76
DF200 (200 ml)	236	4	59	236	4	59	295	5	59

TWO Year Parts Warranty • FIVE Year Vacuum Warranty


Conforms to MDD 93/42/EEC, the Medical Device Directive for the EU. Freezer systems UL/C-UL Listed.

* Without inventory


**Minimum width required for vessel to pass through opening. Footprint may vary. Contact Tech Service for detailed drawings.

Rack Layouts

Square Rack Layout (P)

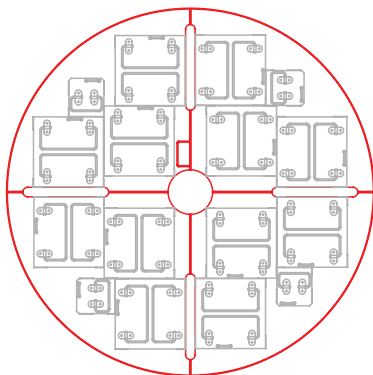

791 OS/U Layout

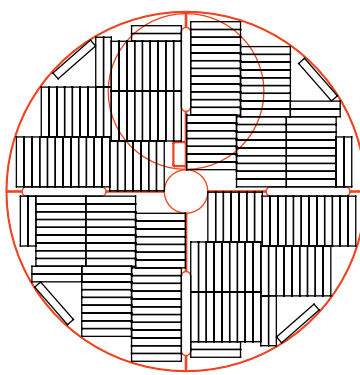

Visit www.chartbiomed.com for more information.

© 2019 Chart Inc. reserves the right to discontinue its products, or change the prices, materials, equipment, quality, descriptions, specifications and/or processes to its products at any time without prior notice and with no further obligation or consequence. All rights not expressly stated herein are reserved by us, as applicable.